



ILLUMINATION INCUBATOR



BSG Series

Shanghai Boxun Medical Biological Instrument Corp
上海博迅医疗生物仪器股份有限公司

Add: No. 599 Zhongqiang Road, Maogang Town,
Songjiang District Shanghai, China

Tel: 021-56980111 Fax: 021-56303876

E-mail: export@shbxyl.com


web: www.sh-bx.com



www.sh-bx.com

ILLUMINATION INCUBATOR >>

BSG Series

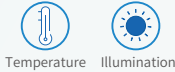

Product Feature

- The shell is made of cold-rolling steel with static spray plastics, the working chamber is made of mirror stainless steel, the shelves can be adjusted optionally;
- PMMA II control system, LCD screen displays temperature, time and illuminance, over-temperature alarm function;
- Observation glass inner door;
- Test hole is equipped on the top, internal diameter is 30mm;
- Well-known international brand compressor, non-CFC refrigerant ;
- Equipped with mechanical lock;
- Integrated type refrigeration system;
- Programmable design, it can be set 30 segments 99 cycles;
- Illuminance can be adjusted in six steps, three-dimensional lighting type;
- Power supply recovery function, data will not be lost when power failure and system halted.



Compressor



Glass Lock



Illumination Board

Model	BSG-250	BSG-300	BSG-400	BSG-800
Volume [L]	250	300	400	800
Temp. Range without Lighting [°C]	0~60	0~60	0~60	0~60
Temp. Range with Lighting [°C]	10~60	10~60	10~60	10~60
Temp. Resolution [°C]	0.1	0.1	0.1	0.1
Temp. Fluctuation at 37°C [°C]	±0.5	±0.5	±0.5	±0.5
Temp. Uniformity at 37°C [°C]	±1.5	±1.5	±1.5	±2
Illuminance [Lux]	0~15000	0~15000	0~15000	0~15000
Power [W]	1300	1500	1500	3300
Timer Range [min]	0~9999	0~9999	0~9999	0~9999
Internal Dimensions W×D×H [mm]	480×480×1100	540×520×1100	554×610×1148	1220×585×1123
Housing Dimensions W×D×H [mm]	710×775×1770	770×815×1770	783×905×1818	1633×910×1820
Number of Shelves (std./max.) [pc]	4/7	4/7	4/7	6/12
Net Weight [kg]	179	203	230	410
Gross Weight [kg]	206	236	250	450
Package Size [mm]	800*910*1960	920*1000*1950	930*1090*2000	1730*1040*2020